

FREEWAY Falls Recovery Sling



The Freeway Falls Recovery Sling allows carers to safely transfer a person that needs assistance immediately after a fall. With a versatile design for both adults and children, the sling has a maximum weight capacity of 200kg.

The Freeway Falls Recovery Sling allows for the fallen person to be rolled onto the sling whilst in a lying position. Once moved onto to the sling, the person can be easily slid to a more open area if required, for the suitable positioning of a hoist. The person can then be hoisted up off the floor and safely into a chair or bed.

It is important for the fallen person to be assessed before any sling application to ensure that injuries haven't be sustained and medical attention isn't required.

PRISM MEDICAL UK


Key benefits and features

- Allows for safe transfer of a person that has fallen on the floor to a chair or bed
- High maximum weight capacity of 200kg (31.5 stone)
- Strong, yet comfortable build
- Clinically proven spacer material for maximum comfort
- Low friction material for comfortable transfer and sliding of person (the yellow side)
- Easy application and removal
- Space saving storage and carry bag (integral)
- Machine washable



How to use the Freeway Falls Recovery Sling

The Freeway Falls Recovery Sling is easy to apply for a qualified carer and the storyboard below illustrates the process in which a person that has fallen is hoisted to a chair.

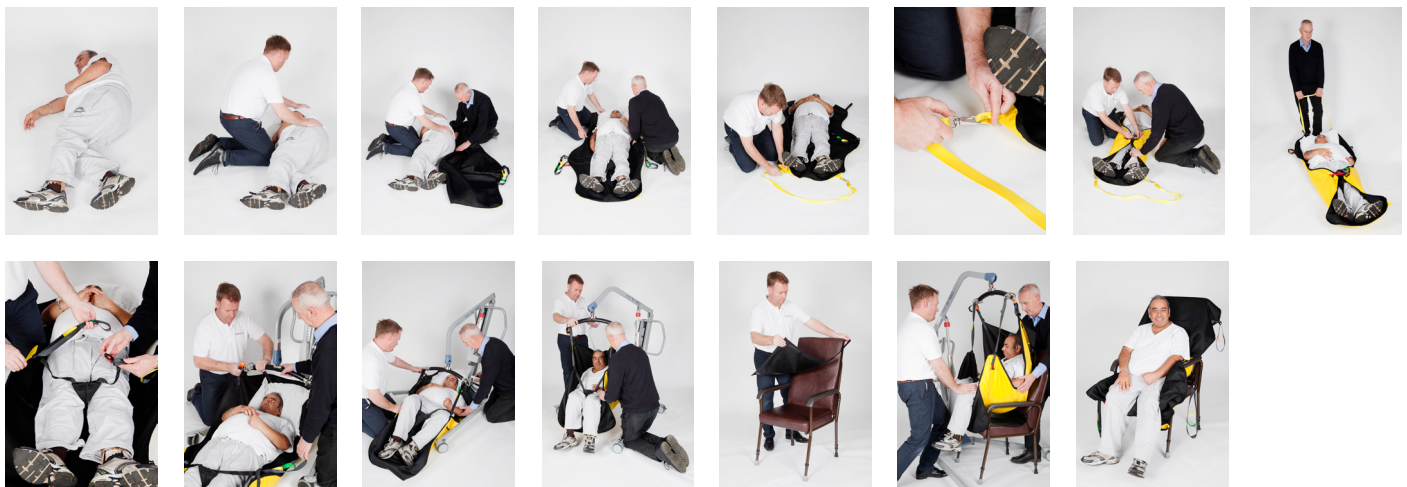

Please refer to the Freeway Falls Recovery Sling User Guide for full instructions on how to safely use this product. The carers must be assessed to be a competent person, trained in the suitable application and fitting of slings.

Contact us:

0845 409 2000 or visit www.prismmedical.co.uk

Registered Office: Prism Medical UK, Unit 4, Jubilee Business Park, Jubilee Way, Grange Moor, Wakefield WF4 4TD
Sling Manufacturing: Prism Medical UK, Millers Avenue, Brynmenyn Industrial Estate, Bridgend CF32 9TD

